

HRSV Post-Fusion glycoprotein F0 Specific ELISA Kit (For Vaccine Development)

Pack Size: 96 tests

Catalog Number: RAS-A204

IMPORTANT: Please carefully read this manual before performing your experiment.

For Research Use Only. Not For Use in Diagnostic or Therapeutic Procedures



INTENDED USE

It was developed for the specific quantitative detection of HRSV Post-Fusion glycoprotein in samples. It is intended for research use only (RUO).

PRINCIPLE OF THE ASSAY

Most in vitro RSV neutralizing antibodies in human sera are directed against the prefusion conformation, but due to its instability the pre-fusion conformation has a propensity to prematurely refold into the stable post-fusion conformation, both in solution and on the surface of the virions. An RSV F protein that has both high expression levels and maintains a stable pre-fusion conformation would, therefore, be a promising subunit vaccine candidate against RSV. To facilitate the RSV-related research, drug trials and vaccine development, a rapid and effective assay kit detecting the specific levels of HRSV Post-Fusion glycoprotein F0 is urgently needed to accelerate the development of RSV vaccines.

This assay kit is used to measure the levels of HRSV Post-Fusion glycoprotein by employing a standard sandwich-ELISA format. The microplate in the kit has been pre-coated with Anti-Post-Fusion glycoprotein F0 (RSV) Antibody. First add the standard samples provided in kit and your samples to the plate, incubate and wash the wells. Then add the HRP-Anti-RSV-F0 Antibody to the plate, incubate and wash the wells. Lastly load the substrate into the wells and monitor color development in proportion with the amount of HRSV Post-Fusion glycoprotein present. The reaction is stopped by the addition of a stop solution and the intensity of the absorbance can be measured at 450 nm and 630 nm. The OD Value reflects the amount of HRSV Post-Fusion glycoprotein bound.

MATERIALS PROVIDED

TABLE 1. MATERIALS PROVIDED

Catalog	Components	Size (96	Format	Storage		
	•	tests)		Unopened	Opened	
RAS204-C01	Pre-coated Anti-Post-Fusion glycoprotein F0 (RSV) Antibody Microplate	1 plate	Solid	2-8°C	2-8°C	
RAS204-C02	Post-Fusion glycoprotein F0 (RSV) Standard	20 μg	Powder	2-8°C	-70°C	
RAS204-C03	HRP-Anti-RSV-F0 Antibody		Powder	2-8°C, avoid light	-70°C, avoid light	

2 / 10

US and Canada:

Tel: +1 800-810-0816

Web: http://www.acrobiosystems.com
E-mail: order@acrobiosystems.com

Asia and Pacific:

Tel: +86 400-682-2521



RA204-EN.01

RAS204-C04	10×Washing Buffer	50 mL	Liquid	2-8°C	2-8°C
RAS204-C05	2×Dilution Buffer	50 mL	Liquid	2-8°C	2-8°C
RAS204-C06	Substrate Solution	12 mL	Liquid	2-8°C, avoid light	2-8°C, avoid light
RAS204-C07	Stop Solution	7 mL	Liquid	2-8°C	2-8°C

REAGENTS/EQUIPMENT NEEDED BUT NOT SUPPLIED

Single or dual wavelength microplate reader with 450 nm and 630 nm filter;

Centrifuge;

37°C Incubator;

10 μL, 200 μL and 1000 μL precision pipettes;

 $10 \mu L$, $200 \mu L$ and $1000 \mu L$ pipette tips;

Multichannel pipettes;

Tubes;

Graduated cylinder to prepare Wash Solution;

Deionized or distilled water to dilute 10×Washing Buffer;

STORAGE

Unopened kit should be stored at 2°C -8°C upon receiving.

The opened kit should be stored per TABLE 1. The shelf life is 30 days from the date of opening.

Note: a. Do not use reagents past their expiration date.

b. Find the expiration date on the outside packaging.

REAGENT PREPARATION

- 1. Bring all reagents and samples to room temperature (20°C-25°C) before use. If crystals have formed in buffer solution, place the sample in a 37 °C incubator until the crystals have completely dissolved and bring the solution back to room temperature before use.
- 2. Reconstitute the provided lyophilized materials to stock solutions with distilled, sterile water as recommended in Table 2 and place the materials for 15 to 30 minutes at room temperature with occasional gentle mixing. Avoid

3 / 10

US and Canada:

Tel: +1 800-810-0816

Web: http://www.acrobiosystems.com

Asia and Pacific: Tel

Tel: +86 400-682-2521

E-mail: order@acrobiosystems.com





vigorous shaking. The reconstituted stock solutions should be stored at -70°C. It is recommended not to freeze-thaw more than 1 times, the packing specification shall not be less than 5 μg.

TABLE 2. RECONSTITUTION METHODS FOR 96 TESTS

ID	Components	Size	Stock Solution Con.	Reconstitution Buffer and
RAS204-C02	Post-Fusion glycoprotein F0 (RSV) Standard	20 μg	100 μg/mL	200 μL water
RAS204-C03	HRP-Anti-RSV-F0 Antibody	20 μg	100 μg/mL	200 μL water

RECOMMENDED SAMPLE PREPARATION

1. Working fluid preparation

1.1 Preparation of 1×Washing Buffer:

Dilute 50 mL 10×Washing Buffer with ultrapure water/deionized water to 500 mL.

1.2 Preparation of 1×Dilution Buffer:

Dilute 50 mL 2×Dilution Buffer with 1×Washing Buffer to 100 mL.

1.3 Preparation of HRP-Anti-RSV-F0 Antibody working fluid:

Dilute HRP-Anti-RSV-F0 Antibody to $0.2~\mu g/mL$ with Dilution Buffer. The prepared working fluid should avoid light. Please prepare it for one-time use only.

2. Preparation of Standard curve

Make serial dilutions of the post-Fusion glycoprotein F0 (RSV) as a Standard curve with Dilution Buffer as recommended in Figure 1.

4 / 10

US and Canada:

Tel: +1 800-810-0816

E-mail: order@acrobiosystems.com



300 µL



Post-Fusion Tubes/ glycoprotein F0 Solution (RSV) Standard Std.-0 Std.-1 Std.-2 Std.-3 Std.-4 Std.-5 Std.-6 Code stock solution 300 μL 10 µL 300 µL 300 uL 300 µL 300 µL 10µL Operating 2000 1.563 0.781 25 12.5 6.25 3.125 Solution 100µg/mL ng/mL ng/mL ng/mL ng/mL ng/mL ng/mL Con. ng/mL

FIGURE 1. PREPARATION OF 1:1 SERIAL DILUTIONS OF THE Post-Fusion glycoprotein F0 (RSV)

6. Add Samples

Dilution

Buffer Vol.

Add 100 μL serially diluted **post-Fusion glycoprotein F0 (RSV)** Standard curve and samples to each well. For blank Control wells, please add 100μL 1×Dilution Buffer. Seal the plate with microplate sealing film and incubate at 37°C for 1.0 hour.

300 µL

790 µL

490 µL

300 μL

300 µL

300 µL

7. Washing

Remove the remaining solution by aspiration, add 300 µL of 1×Washing Buffer to each well, soak for 30 s, remove any remaining 1×Washing Buffer: by aspirating or decanting, invert the plate and blot it against paper towels. Repeat the wash step above for three times.

8. Add HRP-Anti-RSV-F0 Antibody

For all wells, add 100 μL **HRP-Anti-RSV-F0 Antibody (dilute to 0.2 μg/mL)** working solution. Seal the plate with microplate sealing film and incubate at 37°C for 1.0 hour.

9. Washing

Repeat step 4.

10. Substrate Reaction

Add 100 µL Substrate Solution to each well. Seal the plate with microplate sealing film and incubate at 37°C for 20

5 / 10

 US and Canada:
 Tel: +1 800-810-0816
 Web: http://www.acrobiosystems.com

 Asia and Pacific:
 Tel: +86 400-682-2521
 E-mail: order@acrobiosystems.com



ACCO*

min, avoid light.

11. Termination

Add 50 µL Stop Solution to each well and tap the plate gently to allow thorough mixing.

Note: the color in the wells should change from blue to yellow.

12. Data Recording

Read the absorbance at 450 nm and 630 nm using UV/Vis microplate spectrophotometer within 5 minutes.

Note: To reduce the background noise, subtract the value read at $OD_{450 \text{ nm}}$ with the value read at $OD_{630 \text{ nm}}$.

CALCULATION OF RESULTS

- 1. Normal range of Standard curve: R²≥0.9900, detection range: 0.781-25 ng/mL.
- 2. If the OD value of the sample to be tested is higher than the highest standard, the sample shall be diluted with dilution buffer and assay repeated.
- 3. To calibrate absorbance value obtained by the standard curve, the OD value of the sample to be measured is subtracted to the OD value of the blank control. The standard curve is plotted with the standard concentration as x-axis and the calibrated absorbance value as y-axis. Four parameters logistic are used to draw the standard curve and calculate the sample concentration.

PRECAUTIONS

- 1. This kit is for research use only and is not for use in diagnostic or therapeutic procedures.
- 2. The kit should be used according to the instructions.
- 3. Do not mix reagents from different lots.
- 4. Bring all reagents and samples to room temperature (20°C-25°C) before use. If crystals have formed in
- 5. buffer solution, warm to room temperature until the crystals have completely dissolved.
- 6. The kit should be stored at 2°C to 8°C.


TYPICAL DATA

The following data is for reference only. The sample concentration was calculated based on the results of the standard curve

6 / 10





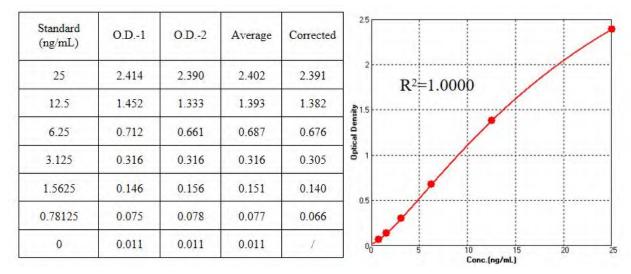

PRECISION

1. Intra-assay Precision

Three samples of known concentration were tested ten times on one plate to assess intra-assay precision.

2. Inter-assay Precision

Three samples of known concentration were tested in three separate assays to assess inter-assay precision.

	I	ntra-assay Precision	n]	Inter-assay Precision	n
Sample	1	2	3	1	2	3
n	10	10	10	3	3	3
Mean (ng/mL)	18.439	4.666	1.834	18.112	4.560	1.812
SD	1.037	0.217	0.107	0.566	0.184	0.038
CV (%)	5.6	4.7	5.8	3.1	4.0	2.1

Note: The example data is for reference only.

RECOVERY

Three samples with different concentrations were tested to calculate the recovery rate.

Sample(n=5)	Average Recovery %	Range %
High	91.0	80.6-105.2
Middle	97.2	88.7-114.2

7 / 10

US and Canada:

Tel: +1 800-810-0816

 $\textbf{Web} \colon \underline{\text{http://www.acrobiosystems.com}}$

Asia and Pacific:

Tel: +86 400-682-2521

E-mail: order@acrobiosystems.com



RA204-EN.01

Low 99.8 94.2-109.3	Low	99.8	94.2-109.3
---------------------	-----	------	------------

LINEARITY

To assess the linearity of the assay, samples spiked with high concentrations were serially diluted with calibrator diluent to produce samples with values within the dynamic range of the assay.

		Cell culture medium (DMEM)	Cell culture medium (1640)	
1:2	Average Recovery (%)	97.2	102.5	
1:2	Range (%)	92.5-105.5	98.0-106.6	
1:4	Average Recovery (%)	93.9	100.3	
1:4	Range (%)	89.2-97.5	92.7-106.9	
1:8	Average Recovery (%)	99.2	107.5	
1:8	Range (%)	92.0-104.9	105.4-109.7	
1.16	Average Recovery (%)	98.0	104.6	
1:16	Range (%)	92.6-111.4	96.0-114.2	

Note: The example data is for reference only.





PLATE LAYOUT

	1	2	3	4	5	6	7	8	9	10	11	12
Α	Std1	Std1							(()	()	
В	Std2	Std2	((\cdots)	()		()	($$ $)$	()	$\left(\right)$	$\left(\begin{array}{c} \cdots \end{array} \right)$	$\left(\cdots \right)$
С	Std3	Std3	()	($$ $)$	(\mathcal{A}	()	()	()	$\left(\cdots \right)$	(($$ $)$
D	Std4	Std4	(($$ $)$	(()	($)$	()	()	($ $	($)$
E	Std5	Std5	<u></u>	($)$	\(\tag{\tag{\tag{\tag{\tag{\tag{\tag{		$\langle \cdots \rangle$	$\langle \cdots \rangle$	\\	()	$\langle \cdots \rangle$	()
F	Std6	Std6		($)$	\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\	$\stackrel{\dots}{\searrow}$	$\stackrel{\dots}{\searrow}$	()		<u></u>	()	()
G	Blank	Blank	\;;;;	$\langle \cdots \rangle$	\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\	$\stackrel{\dots}{\nearrow}$	$\langle \rangle$	\ <u>\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\</u>	\\	\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\	\ <u>\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\\</u>	()
Н	Blank	Blank	((-)	· · ·)		(()	(()	()	()

Note: Blank is a Blank Dilution Buffer hole.

TROUBLESHOOTING GUIDE

Problem	Cause	Solution		
Poor standard curve	* Inaccurate pipetting	* Check pipettes		
Lauga CV	* Inaccurate pipetting	* Check pipettes		
Large CV	* Air bubbles in wells	* Remove bubbles in wells		
Walle day of	* Plate is insufficiently washed	* Review the manual for proper wash.		
High background	* Contaminated wash buffer	* Make fresh wash buffer		
Very low readings across the	* Incorrect wavelengths	* Check filters/reader		
plate	* Insufficient development time	* Increase development time		



RA204-EN.01

Samples are reading too high, but standard curve looks fine	* Samples contain cytokine levels above assay range	* Dilute samples and run again			
		* Assay set-up should be continuous - have all standards			
	* Interrupted assay set-up	and samples prepared appropriately before commencement of theassay			
Drift	* Reagents not at room temperature	* Ensure that all reagents are at room temperature before			
		pipetting into the wells unless otherwise instructed in the			
		antibody inserts			